# S · E · Y · S · S · O

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
## SET ELEMENTS

- 1 irrigator handle
- 1 water tank (150 ml)
- 3 nozzles
- 1 inductive charging station

## DESCRIPTION OF SEYSSO CARBON TRAVEL IRRIGATOR

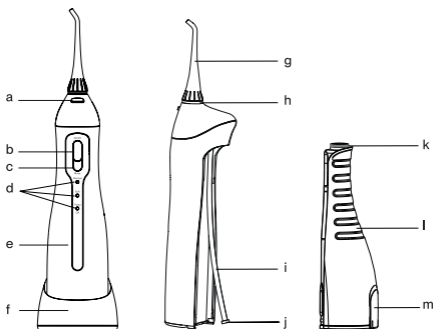- |                          |                        |
|--------------------------|------------------------|
| a. Nozzle release button | g. Nozzle              |
| b. On/Off button         | h. Identification Ring |
| c. Mode button           | i. Suction Tube        |
| d. Mode indicator        | j. Filter              |
| e. Handle                | k. Water Tank Gasket   |
| f. Charging station      | l. Water Tank          |
|                          | m. Water Tank Cap      |

### NOTE:

There may be some liquid remaining inside the main appliance before use (tank, water suction tube, etc.), this is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.

## PLUGGING AND CHARGING

### CHARGING THE DEVICE (FIG. 1)

- Plug the charger into a power supply.
- Put the device on the charger in the right position.
- The CH indicator keeps blinking while charging and turns off as soon as the device is fully charged.

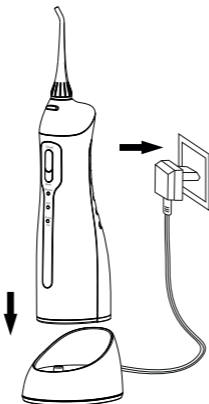
Fig.1

### NOTE:

If the CH indicator keeps blinking during or after use, that means the battery needs to be recharged. Full recharging takes approximately 24 hours.

- Before the first use and after the irrigator has not been used for at least 3 months, charge the device fully for 24 hours. The fully charged device can operate for up to 70 minutes. This provides 14 days of use without recharging if the irrigator is used twice a day for 2.5 minutes.
- Charge the device only when indicated. Additional charging may reduce battery life. Discharging the battery fully is not recommended.

## USAGE OF THE ORAL IRRIGATOR SEYSSO TRAVEL

- Attach the nozzle to the main unit. (Fig.2)
- Open the water tank cap to fill the tank with water, then close the cap. The temperature of the water should not exceed 40°C. (Fig. 3)
- Alternatively, slide the water tank down off the main unit, fill the tank with water not warmer than 40°C, then reattach the tank to the main unit. (Fig. 4)
- Press the Mode button to select the mode you want to use. You can change the working mode while using the irrigator (Fig. 5)

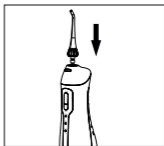
Fig. 2

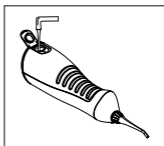
Fig. 3

**NORMAL:** for normal cleaning (6,2 bar).

**SOFT:** for gentle cleaning of sensitive gums (2,7 bar).

**PULSE:** for cleaning and massaging the gums (2,7-6,2 bar).

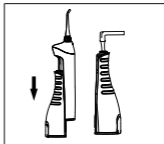
Fig. 4

**WARNING:** Use Carbon Pocket Nozzle on soft mode. High pressure can damage it.

- Hold the unit in an upright position with the nozzle pointing toward your teeth.
- Keep your mouth slightly open so that the water can flow out.
- Press the „On/Off” button to start irrigating. (Fig. 6)

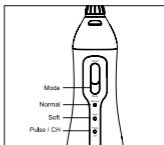
Fig. 5

- ① Direct the water jet flow at a 90° angle to the teeth and gums.
- ② Move the water jet slowly along the teeth.
- ③ Move the water jet along the gum line with the nozzle at a 90° angle to the gum tissue.
- ④ Clean between the teeth and around any braces, crowns or bridgework thoroughly.

- Press the „On/Off” button and then press the nozzle release button to replace the nozzle. (Fig. 7)
- Once the irrigation has been completed, press the „On/Off” button again to turn off the device.

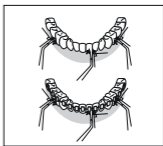
Fig. 6

#### AFTER USING THE APPLIANCE (FIG. 8)

- Open the water tank cap and dispose of the water in the tank.
- Press „On/Off” button to turn on the unit, and drain the water inside the appliance.
- Press „On/Off” button again to turn off the unit, and close the water tank cap.
- Wipe the device with a dry cloth, a soft cloth.

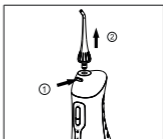
Fig. 7

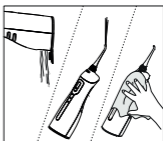
Fig. 8

#### NOTE:

To prevent bacterial growth, make sure there is no remaining water inside the unit. If you are not going to use this appliance for a long period, be sure to wipe and dry the appliance before storing it.

## COMPATIBILITY OF NOZZLES

The SEYSSO Carbon Travel irrigator is compatible with SEYSSO Oxygen/Carbon Standard, SEYSSO Oxygen/Carbon Pocket and SEYSSO Oxygen/Carbon Orthodontic nozzles. We do not guarantee the correct operation of the device in case of using other nozzles than the ones mentioned above.

## FUNCTIONS OF THE NOZZLES

### ORTHODONTIC NOZZLE

The SEYSSO Carbon Orthodontic nozzle is designed for people with fixed braces.

- Place the SEYSSO Carbon Orthodontic nozzle on the handle.
- Set a low water pressure (Soft Mode), and then gradually increase it as needed by selecting Normal or Pulse mode.
- Direct the water stream at a 90° angle to your teeth, at the gum line. Start irrigating from the last molar.
- Slowly glide the nozzle along the gumline, stopping briefly to clean interdental spaces, and areas around braces more thoroughly.
- Continue irrigating until the inside and outside of all teeth are precisely cleaned.

### POCKET NOZZLE

The SEYSSO Carbon Pocket nozzle is specifically designed to clean the gingival pockets. To use the Carbon Pocket nozzle, ALWAYS set the device to the lowest pressure value (Soft Mode).

- Place the nozzle against a tooth at a 45° angle with the soft, rubber tip, facing the gingival pocket.
- Turn on the device, and move the nozzle along the gumline until all gum pockets are cleaned.

**WARNINGS! READ IT BEFORE USING THE APPLIANCE!**

### SAVE THIS MANUAL FOR FURTHER REFERENCE

- For domestic use only.
- The irrigator is equipped with a waterproof handle. It is safe and designed to be used in the bathroom. The appliance has a IPX7 protection rating. This means that the irrigator is resistant to accidental submersion in water up to the depth of 1 metre for up to 30 minutes.

- Remember: the irrigator is not suitable for operating underwater.
- The product is not intended for repair on your own.
- Store the charger in a dry and safe place.
- The appliance is only to be used with the power supply unit provided with the appliance.
- Straighten the cable before using it.
- In the event of a failure, do not dismantle the appliance or replace any parts, as there is a risk of being electrocuted. Contact the service centre SEYSSO.
- A damaged power cable must be replaced with a new one. Contact the service centre SEYSSO.
- If dropped, replace the nozzle before its next use. Even if the damage is not visible.
- The device should not be disassembled except prior to disposal, in which case the battery should be removed.
- The battery is to be disposed of safely.
- Refer to "Disposal" on page 23 for more information about removing the battery.
- The appliance must be disconnected from the supply mains when removing the battery.
- Do not fill the tank with water warmer than 40°C.
- Do not use the device if you lost the nozzle.
- Keep the nozzle away from small children, as it is small enough to be swallowed and can be lodged in the throat.
- Do not fill the device with mouthwash as it may damage the equipment by cracking the water tank.
- If during or after irrigation you feel pain or discomfort, or if you are concerned that the device may have any negative impact on your health, please consult a doctor.
- If in the last two months you have undergone any surgery of the oral cavity, please consult a doctor before using the device.
- Do not press the nozzle too firmly against your teeth or gums. Doing so may injure your gums.
- The irrigator is designed to clean the interdental spaces, gums and tongue. Do not use it in any other way.
- Use the SEYSSO Carbon Pocket nozzle only on the Soft Mode. High pressure can damage the tip.

## MAINTENANCE AND CLEANING

- Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device.
- Do not wash using hot water of over 50°C.
- Keep the device away from high temperature areas and direct sunlight.
- The irrigator should be descaled once a month using a cleaning agent designed for this purpose.

### MAIN UNIT

- Remove the dirt by wiping using a cloth with liquid soap and wipe off using a dry cloth, a soft cloth.
- Do not submerge the main unit in water to clean.
- Dry after washing with water.

### INDUCTION CHARGER

Check the charger plug and cable at regular intervals, wipe the adapter plug with a dry cloth, a soft cloth.

### NOZZLE

- Wash with water and wipe with a dry cloth, a soft cloth.
- Replace the nozzle every 6 months. Replace the Orthodontic Nozzle and the Pocket Nozzle every 3 months.

### WATER SUCTION TUBE

- Wash with water and wipe with a dry cloth, a soft cloth.
- Do not bend, pull or twist the water suction tube.

### WATER TANK

- Wash with water and wipe with a dry cloth, a soft cloth.
- Dry the tank when not using for a long period.

**DANGER WARNING**

In order to reduce the risk of being electrocuted:

- Always disconnect the appliance from the power source after use.
- Do not use during bathing.
- For more detailed instructions related to storage and use of the appliance, as well as the attached accessories, go to SEYSSO.

In order to reduce the risk of catching fire, short circuit or injury:

- Disconnect the charger from the power source when charging is complete.
- The product is not a toy. Take special caution when the irrigator is being used by a child or in the presence of a child.
- Use it in accordance with its intended use.  
Use only accessories recommended by the manufacturer.
- Do not use the charger if the power cable or plug is faulty.  
If the waterflosser does not work properly after being dropped, left in the open space, or submerged in water, send it back to the place of purchase.
- If the charger falls into water, do not try to take it out, instead unplug it immediately.
- Do not stretch the charger power cable, jam it with a door, wrap it around sharp edges. Do not switch on the irrigator while it is charging. Keep the appliance components away from hot surfaces.
- Do not disconnect the charger from the power source by pulling on the cable. Use the plug.
- Do not keep the plug with wet hands.
- Never use the product if you have problems concentrating and if you experience drowsiness.
- Do not use the appliance outdoors or in the places where it can be exposed to aerosols or pure oxygen.
- Plug the appliance into the socket. Read its specifications.

It is not recommended that children under the age of 8 years use the device. Children and persons with limited sensory, motor or mental abilities should use the appliance under the supervision of a guardian, or should be properly trained in the scope of correct and safe use of the appliance.

## ADDITIONAL INFORMATION

- Do not share your nozzle with other users.
- Do not wind the power cord of the charger too tightly.
- Do not put any metal object, e.g. coins, clips, on the charger.
- Unplug the charger when you know it will not be used in the immediate future.
- The gap near the reservoir mount is intentional, and does not have influence on the tightness of the handle.

## DISPOSAL

The rechargeable battery inside your appliance is easy to remove for recycling. At the end of the product's life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below. During the following operation, be sure to protect your eyes, hands, fingers, and the surface you work on.

- Remove 2 rubbers feet on the bottom of the handle. (① on Fig. 9)
- Insert a screwdriver into the slot located at the bottom of the handle and turn counter-clockwise to release the bottom cap. (② on Fig. 9)
- Slide the water tank down off the device and unscrew the two screws next to the suction tube. (Fig. 10)
- Open the top cover and take the internal part out of the unit, then remove all of its screws. (Fig. 11)
- Cut off the battery line and pull the batteries out of the internal part. (Fig. 12).

If you have recycling questions, please contact your local waste management office.

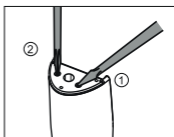
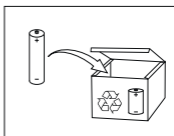
Fig. 9

Fig. 10

Fig. 11

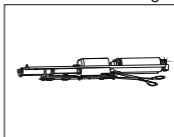
Fig. 12

## TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
It does not work.	You've just purchased the product, or it has been left unused for more than 3 months.	Charge the unit for about 24 hours.
It only works for a few minutes, even after charging.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the unit for about 24 hours.
Water pressure is very low.	There is dirt stuck in the nozzle. Hard water can cause the nozzle to clog.	Replace the nozzle.
	The nozzle is deformed.	
Water does not come out.	The tank is empty.	Fill up the tank.
	The appliance is tilted too much while in use.	Use with the appliance in a vertical position.

## SEYSSO CARBON TRAVEL ORAL IRRIGATOR SPECIFICATION

NAZWA PRODUKTU	SEYSSO Carbon Travel Oral Irrigator
Voltage	100-240 V~, 50/60 Hz
Power	0,7-1,7W
On/Off Switch	Button
Display	LED
Displayed information	Mode + Charging Status Indicator
Maximum Water Pressure	6,2 bar
Water Tank Volume	150 ml
Usage Time	Up to 70 minutes
Charging	Inductive
Charging time	Approx. 24 hours
Battery	Li-Ion
Discharged Battery Indicator	Flashing diode
Water Tightness	IPX7

## 2-YEAR LIMITED WARRANTY

SEYSSO gives the warranty for the Waterflosser for a period of two years, from the date of purchase, against defects due to faulty workmanship or materials. Remember to keep your original purchase receipt for the duration of the warranty period to validate your warranty. Within the valid limited warranty period of this device, the Customer may contact SEYSSO service consultant only when the Product purchased becomes defective when properly used. Our service will replace the device free of charge.

### THE WARRANTY DOES NOT COVER:

- ① Cosmetic damages.
- ② Damage due to:
  - a) accident, misuse, negligence, commercial use or modifications of the product;
  - b) improper operations or maintenance of the Product;
  - c) connection to improper voltage supply;
  - d) attempted or unauthorized repair made by party other than SEYSSO.

Any attempt to open or dismount the device (or its accessories) will void the warranty. Claims under the warranty must be supported with reasonable evidence. Claim date must be within the warranty period. Shipping costs are non-refundable. Do not under any circumstances, scratch or scrape-off the handle coating. All permanent scratches will be considered cosmetic damage and are not covered by warranty. The warranty covers the parts that affect the Product operation.

The instructions apply to models with the numbers:

SC.WF.010.BLK.1/POI01

**MANUFACTURER:**

WM BRANDS SP. Z O. O.

UL. DZIUPIŁI 12

02-454 WARSZAWA, POLSKA

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**CONTACT:**

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www.seysso.pl

Charger input voltage: 100-240 V~, 50/60 Hz

Handle charging current: 60-95mA

Water resistance (body, charger): IPX7

Protection against the effects of immersion in water at a depth of up to 1 meter for 30 minutes Irrigators have a built-in lithium-ion battery.

The product meets the requirements of the European Union (EU) "New Approach" directives.

The packaging and instruction manual are recyclable.

This device is marked in accordance with the Polish Act on Waste Electrical and Electronic Equipment and the WEEE Directive 2012/19/EU with the symbol of a crossed-out waste bin. Such marking informs that this equipment, after its use, cannot be disposed of together with other household waste. The user is obliged to hand it over to entities collecting waste electrical and electronic equipment. Collectors, including local collection points, stores, and municipal units, create an appropriate system allowing for the disposal of this equipment. Proper handling of waste electrical and electronic equipment helps to avoid harmful consequences for human health and the natural environment, resulting from the presence of hazardous components and improper storage and processing of such equipment.

Symbol of Class II protection designation for electrical equipment.

For domestic use only.

S·E·Y·S·S·O

REGISTER YOUR  
PRODUCT  
STREFASEYSSO.PL

A series of several thin, overlapping circles of varying diameters, centered around the text 'Smile matters', creating a sense of motion or a ripple effect.

**Smile matters**

[WWW.SEYSSO.PL](http://WWW.SEYSSO.PL)